

HeartX[®] Recorder

Holter ECG-Recorder



Portable and wireless ECG (no lead wire)

Precise measurements with 3 channels

Waterproof (IP67)

Only 3 min charging for 24h measurement cycle

Up to 7 days continuous recording without charging

Convenient charging and downloading via 2-in1 docking station

Simplified Holter Monitoring for Everyday Use

The HeartX Recorder is designed for ease of use and offers a comfortable and convenient solution for both patients and clinical staff. With its lightweight, wireless design, the device ensures continuous ECG monitoring without interfering with daily activities. This makes the HeartX Recorder ideal for long-term use. The device can be used with standard electrodes or special patches to ensure optimal data quality for patients with different indications.

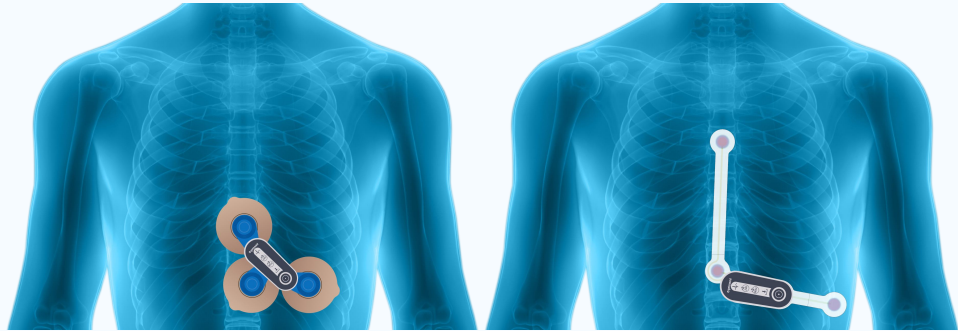
The simple workflow is characterised by fast downloading of the recording and simultaneous charging of the device via the docking station. The magnetic connection of the docking station ensures that the recorder can be placed effortlessly and securely.

HeartX[®] Recorder

Technical Data

Weight	35 g
Dimensions	80.95 x 34 x 14.6 mm [3.19 in x 1.34 in x 0.57 in]
ECG Channels	3
Recording Duration	7 days
Battery Lifetime	5 years
Charging Time	3 minutes for a 24 hour recording
Degree of Ingress Protection	IP 67
Pacemaker Pulse Detection	Yes
Lead-Off Detection	Yes
Operating Temperature	5 °C to 45 °C

Standard Elektrodes and Patches

Different Hookup Options	
Biocompatibility (ISO 10993)	Yes

Figures in original size



GETEMED Medizin- und Informationstechnik AG
Oderstr. 77 / 14513 Teltow / Germany
Telephone +49 3328 3942-0
Telefax +49 3328 3942-99
info@getemed.de / www.getemed.de
Note: Technical data subject to change.